

- Legend:**
- ☑ Good, Do
 - ☹ Don't do

Autoclave Safety

For Liquids:


☑ **ONLY** Pyrex™ or Kimax™



☑ Foam Stopper

☑ Foil Cap


Autoclavable lid **very** loosely threaded ☑

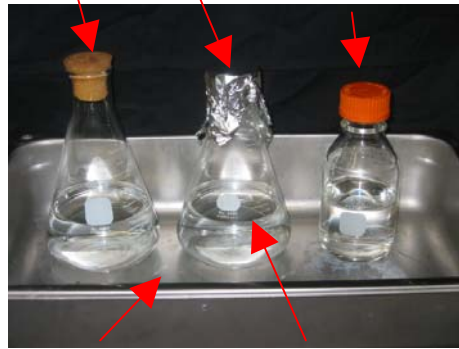

☑ metal tray ☑ volume <50%

↓ ↓ **!!!! DANGER !!!!** ↓ ↓

☹ cap will melt and act like a stopper
 ☹ glass **NOT** borosilicate:
 (likely to crack)

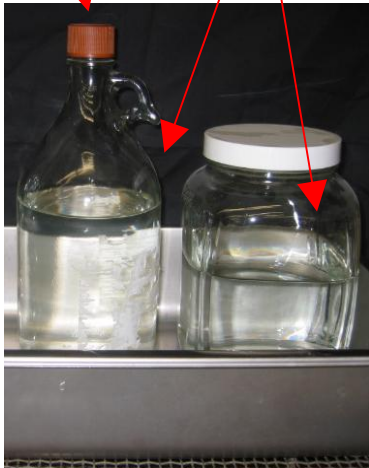

☹ rubber stopper:
 (flask can explode)

☹ too full
 (will boil over)



☹ no metal tray
 (can spill on you)

Solid Biohazardous Waste

- ☑ clear autoclave bag; top open
- ☑ biohazard label
- ☑ metal tray



☹ biohazard bag ☹ top sealed
 (can't go in garbage) (poor heating)



☹ green garbage bag ☹ no metal tray
 (will melt) (big mess)

